


The SixFix® Circular Fixation System

is a platform that can address extremity applications in trauma, foot and ankle, and deformity correction for both adults and pediatrics. The SixFix® System introduces a unique design in which one compact pre-assembled unit of six (6) struts are ready for sterilization and use right out of the box.

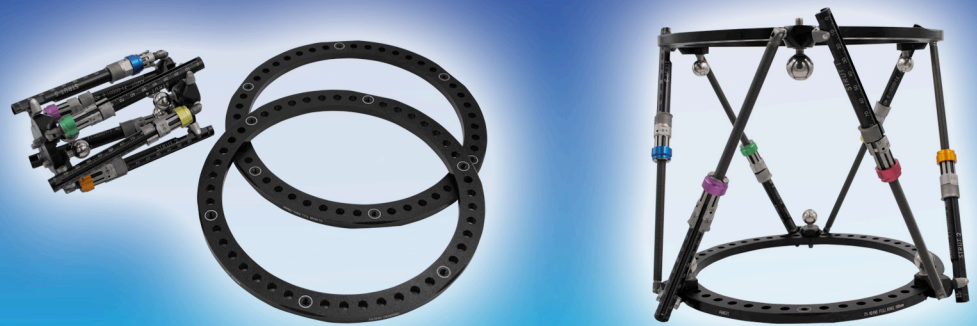

"Each strut in the pre-assembled SixFix® strut assembly is able to extend from 100mm to 300mm. Because I did not have to be concerned about switching out struts intra-operatively as I have with other designs, the application of the frame was fast and easy. At the end of the surgery, I was able to maintain the slim profile of the strut assembly by removing those strut add-ons not required to achieve the final correction."

Stephen J. Frania, DPM, FACFAS

Catalog Information

- 71-00001** Strut Assembly, 100mm - Base Unit
- 73-00003** Strut Assembly, 50mm - Base Unit
- 72-00001** Strut Assembly, 100mm - Base with Add-ons
- 73-00001** Strut Assembly, 50mm - Base with Add-ons
- 71-00002** Strut Screw Add-on, 100mm
- 73-00002** Strut Screw Add-on, 50mm



- 71-00010** Reduction Strut, Short
- 71-00011** Reduction Strut, Medium

Rings

Full Rings

- 71-10075** Full Ring, 75mm
- 71-10095** Full Ring, 95mm
- 71-10115** Full Ring, 115mm
- 71-10125** Full Ring, 125mm
- 71-10135** Full Ring, 135mm
- 71-10145** Full Ring, 145mm
- 71-10155** Full Ring, 155mm
- 71-10165** Full Ring, 165mm
- 71-10175** Full Ring, 175mm
- 71-10185** Full Ring, 185mm
- 71-10195** Full Ring, 195mm
- 71-10215** Full Ring, 215mm
- 71-10225** Full Ring, 225mm

2/3 Rings

- 71-58095** 2/3 Ring, 95mm
- 71-58115** 2/3 Ring, 115mm
- 71-58125** 2/3 Ring, 125mm
- 71-58135** 2/3 Ring, 135mm
- 71-58145** 2/3 Ring, 145mm
- 71-58155** 2/3 Ring, 155mm
- 71-58165** 2/3 Ring, 165mm
- 71-58175** 2/3 Ring, 175mm
- 71-58185** 2/3 Ring, 185mm
- 71-58195** 2/3 Ring, 195mm
- 71-58215** 2/3 Ring, 215mm
- 71-58225** 2/3 Ring, 225mm

Foot Rings & Arches

- 71-6100S** Foot Ring, 150mm, Short
- 71-6100L** Foot Ring, 150mm, Long
- 71-6120S** Foot Ring, 180mm, Short
- 71-6120L** Foot Ring, 180mm, Long
- 71-90150** Foot Arch, 150mm
- 71-90180** Foot Arch, 180mm
- 71-6130S** Foot Rocker, Short
- 71-6130L** Foot Rocker, Long

Assembly

Inserting a Strut Screw Add-On

1. Remove Shoulder Screw. ▼



2. Insert Strut Screw Add-On into the open end of the Strut Body. ▼



3. Place the 3/32" hex driver into the open end of the Add-on.
4. Tighten until the Torque Limiting Handle clicks once.
5. Replace the Shoulder Screw with the Torque Limiting Handle until it clicks once.

Attaching Rings to a Strut Assembly

1. Select the rings to use.
2. Loosen the Locking Ring to allow the struts to extend freely.
3. Remove the Strut Assembly from the Strut Manager.
4. Hand-tighten the Fiducial Bolts into alternating circled holes on both rings.
5. Final tightening is achieved using the Hex Driver, 5mm.



AMDT Holdings, Inc.
328 Poplar View Lane East
Suite 2
Collierville, TN 38017

info@amdtholdings.com
(901) 853-4366

AMDT™

SixFix® Circular Fixation System and the Deformity Analysis & Correction Software (DACs) are FDA cleared in the US.

© 2021 AMDT Holdings, Inc. 78-00004 REV C - 8/2021